

PENGUIN DUAL FEEDING SYRINGE SET

Product Manual English

Product Name: Penguin Dual Feeding Syringe Set
Materials: Nib Tip, Spoon Tip, Syringe Body - Medical-Grade Silicone
 Measuring Cup, Storage Case - PP
 Tube Cleaning Brush - Nylon, PP, Stainless Steel

Dimensions: 1 x Silicone Nib - 40 x 40 x 53mm
 1 x Silicone Feeding Spoon - 40 x 40 x 58mm
 1 x Silicone Syringe Body - 36 x 35 x 70mm
 1 x PP Measuring Cup - 40 x 40 x 45mm
 1 x PP Storage Box - 136 x 83 x 55mm
 1 x Cleaning Brush - 135 x 21 x 6mm

Volume: Syringe Body: 20ml
 Nib Tip: 0.5ml
 Measuring Cup: 20ml

Recommended Age: 0 month+

Product Description:

Feed your baby precise quantities of liquids, purées or medications with the Penguin Dual Feeding Syringe Set! This combo has three different ways to measure your doses, making it practical as well as adorable. The two soft, medical-grade silicone tips are designed to suit different needs, making it incredibly versatile. The dual-port nib tip is ideal for feeding thinner liquids to younger babies, as the ports direct liquid into their cheeks rather than their throats, preventing choking and spit-ups. The feeding spoon, on the other hand, is perfect for everything from thin liquids right through to thicker purées. The X-cut opening allows you to easily control the flow of food onto the shallow easy-feed spoon. This set also includes a measuring cup to quickly and accurately prepare liquids, and a dedicated cleaning brush to return your syringe to a sparkling clean state!

Features:

- Effortless one-handed operation – simply squeeze to fill and control flow rate.
- Precise measurement marks for accurate dosing.
- Two feeding tips to suit different needs.
- Rounded pacifier-shaped nib tip is more readily accepted by babies.
- Nib tip features dual-angled ports to send liquid into the cheeks to prevent choking and spitting up.
- Nib tip's soft shield ensures comfortable feeding without choking.
- Feeding spoon with X-cut valve easily controls flow and dispenses thinner or thicker textures.
- Flexible, transparent syringe body to track intake.
- Tips & body are made from soft and safe 100% medical-grade silicone.
- 20ml PP measuring cup for precise dosing and keeping syringe hygienic between uses.
- Safe in dishwashers, sterilisers and boiling water - simply disassemble to clean.
- BPA, PVC and phthalate-free.

Cleaning and Care:

Before first use, disassemble and clean the product, then place the silicone components in boiling water for 5 minutes. This is to ensure hygiene. Wash before and after each use. Always sterilise after any illness. Although dishwasher safe, we recommend hand washing with the cleaning brush in warm soapy water, then rinsing and drying thoroughly. To sterilise the syringe, measuring cup & storage box, use a steam steriliser or boil in water for 2-3 minutes. Brush **NOT** suitable for microwaves, sterilisers or boiling water due to nylon and stainless steel components. Do not use any bleach-based cleaners or tablets to clean or sterilise your Haakaa products. Silicone may absorb the colours of liquids it comes in contact with. To minimise discolouring, wash immediately after use. Store product in a cool, dry place and avoid direct sunlight.

PLEASE NOTE: Always follow your doctor's or the manufacturer's dosage recommendations. Check measured dosage at eye level on a flat surface. Solid food is recommended for babies 4-6 months+. Please consult your medical professional if introducing solids earlier.


Instructions:



Sterilise syringe and included measuring cup before use. Select the tip that best suits your needs.




Method 1 - Nib Tip. Prepare liquids in the measuring cup. Draw up liquid by squeezing the frosted area of syringe body, submerging the tip and gently releasing. Check the amount again.

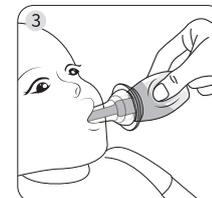

Introduce the tip to the corner of your baby's mouth, and after they start to suck, slowly squeeze the body.



Method 2 - Carefully pour the desired liquid or purées directly into syringe body. Snugly push the tip onto the body, ensuring there are no gaps. Check the amount again.

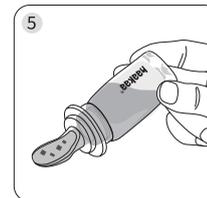

Slowly squeeze until spoon is filled to spoon-feed your baby.



Cleaning: after feeding, disassemble and wash thoroughly with the dedicated brush, according to instructions provided.

NOTE: Please do not forcefully squeeze or shake the body at any time. Test flow strength prior to each use to avoid bubbles forming in the syringe body.

Caution:

For your child's safety and health

WARNING!

- ⊙ Continuous and prolonged sucking of fluids will cause tooth decay.
- ⊙ Always check food temperature before feeding.
- ⊙ Check product condition regularly. Throw away at the first sign of damage or weakness.
- ⊙ Keep components not in use out of the reach of children.
- ⊙ To avoid the risk of strangulation, never attach to cords, ribbons, laces or loose parts of clothing.
- ⊙ Never use tips as a soother.
- ⊙ Always use this product with adult supervision.
- ⊙ Do not store near any sharp objects.
- ⊙ This product is not a toy. Do not use this product for anything other than its intended use.
- ⊙ Inspect tips before each use by pulling them in all directions to ensure no damage exists.
- ⊙ Do not leave this product in direct sunlight or heat, or leave in non-bleach-based sterilising solution for longer than recommended, as this may weaken the product.
- ⊙ Do not fill syringe beyond top measurement mark.
- ⊙ Test flow strength by squeezing on syringe body prior to each use.